

PalmScan™ Mobile USB B-Scan

MMD's Portable and Versatile Ultrasound Imaging

Turn any PC* with a USB port into a high resolution B-Scan



Solutions that Fit


PalmScan™ Mobile USB B-Scan

Turn any PC into a B-Scan device capable of high-resolution ultrasonic imaging. Requires USB 2.0 Interface or higher

The Mobile B-Scan is MMD's portable and extremely versatile ultrasound device. The B-Scan is easy to use with custom settings and a patient database.

The Mobile B-Scan is not only convenient but efficient and accurate with high resolution imaging that gives over 2048 sample points per line. Mobile B-Scan has a probe frequency of 12/15MHz, a 3 to 10 penetration depth and with its dual calipers you can easily measure the structure within the eye.

Printing and capturing images of the eye have never been easier. Simply select desired images for your report and generate reports that can be sent to your existing office printer or store movies or JPEG images right to your PC's hard drive. Adjustable size cine buffers with a maximum of 512 buffers ensure that you always capture your best images.

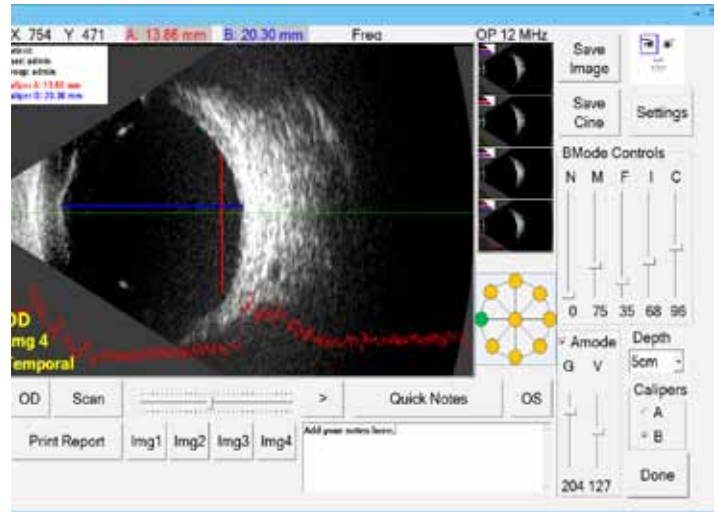

FEATURES

Easy to Use

Custom settings for multiple users.

Patient Database

Virtually unlimited in size depending on PC memory capacity.

High Resolution Imaging

2048 samples per line

Dual Calipers

Easily measure structures within the eye.

Printing Capability

Generate professional pdf reports that can be printed to your existing office printer or transferred to your EMR

Single Frame and Cine Modes

Save movies or JPG images with adjustable size cine buffers.

Zoom Capability

Directly zoom into the structures within the eye

Pulse Frequency 12/15 MHz

3-10 cm penetration depth.

A-Mode Display

Vector selection tools and A-Mode display

Quick Notes

Rapid Note Taking Software

BENEFITS

- Handheld & Mobile
- USB Powered
- Extremely Convenient

SPECIFICATIONS

Depth	3-10 Cm
Probe Frequency	12 - 15 MHz
Samples per line	2048
Weight	5.5 ounces
Power Source	USB



Proudly made in the U.S.A.



A New Generation of Ophthalmic Convenience and Efficiency

Visit Micro Medical Devices at www.micromedinc.com or contact your Authorized Micro Medical Devices Distributor.

IMPORTANT NOTICE

For your own protection, only purchase and service your Micro Medical Devices Inc. devices through an Authorized MMD Distributor.

ISO certified, GMP

Micro Medical Devices, Inc. products are designed and manufactured under quality processes meeting ISO requirements